

Wavelength	Two models, 808nm or 976nm
Power Output	976nm: 0.1-16W@pulse mode, 0.1-10W@CW mode
	808nm: 0.1-12W@pulse mode, 0.1-7W@CW mode
Pulse Duration	10µs-0.9s
Pulse Repetition Rate	Up to 50kHz
Tip Diameter	200μm, 300μm, 400μm
Power Supply Input	100-240VAC, 50-60Hz
Battery Pack	Li-ion Rechargeable, Last for one week with normal use
Size / Weight	$138 \times 142 \times 168$ mm, 1.3 kg

#### **APPLICATIONS**

	Gingivectomy/Gingivoplasty	
	Implant Uncovery	
	Operculectomy	Periodont
	Frenectomy/Frenotomy	
Surgery	Hemostasis	
	Herpes	Endodor
	Aphthous Ulcers	Bio-Stimu
	Fibroma	
	Abscess Excision	Cosme

	Periodontal Pockets
	Peri-implantitis
Periodontology	Gingivitis
	Troughing
	Sulcular Debridement
Endodontics	Pulpotomy
Dr. Cr. Lu	Temporary Pain Relief
Bio-Stimulation	TMJ Disorders
Cosmetics	Tooth Whitening
Cosmetics	Crown Lengthening

# LAZON MEDICAL LASER CO., LTD.

Tel: +86-24-31610402
Fax: +86-24-31610407
Email: service@lazonlaser.com
Web: www.lazonlaser.com



# SOLASE DENTAL DIODE LASER



SOLASE Dental Diode Laser www.lazonlaser.com



#### **LASER**

Less bleeding in the tissue

Painless, no need for anesthesia

Simple and efficient operation. Time-saving, High precision, No noise, No vibration

Operation is safe to metal such as implant.

Little effect on the surrounding tissue

Low probability of cross infection with sterilizing effect

Faster postoperative tissue healing. Small scar with good appearance

With pain relief effect, slight discomfort after

#### **ELECTRIC KNIFE OR SCALPEL**

More bleeding by scalpel

Strong sense of pain by scalpel, need for anesthesia

Difficult to operate. Time-wasting and low precision. Electric knife is dangerous to the device on patient

Electric knife operation is forbidden around impant

Obvious contraction of surrounding tissue.

Bigger heat affected area

Antibiotics need to be taken after operation by scalpel

Slow healing with obvious scar

Severe postoperative swelling with strong discomfort

# TO BE YOUR SMARTEST ASSISTANT

#### **SELF-LEARNING**

The APP includes built-in operating instructions, user manual and instructional video.

#### **PATIENT RECORD**

To save your patient information, treatment settings, photos, video, notes and therapeutic effect evaluation.



#### **APP UPDATING**

Update the APP in the first time to get the latest features for your laser









#### **UNLIMITED CUSTOM PRESET**

The user can save an unlimited number of custom presets.



#### **MULTI-ACCOUNT**

Every dentist can have their own account protected by their own password.



#### CLOUD BACKUP

All data can be backed up to the Cloud and saved in private.



#### **SMART TIP INITIATION**

The APP identifies the indication and automatically turns on or off the tip initiation.



**OFFLINE MODE** 

Check your data anywhere

without the laser by your side.

www.lazonlaser.com



WIDER RANGE OF POWER FROM 0.1-16W EXPANDS **SCOPE OF INDICATIONS** 

> With a power range from 0.1W to 16W, the SOLASE Laser can be used for more than 30 kinds of oral diseases including surgery, periodontics, endodontics, dental implantology, aesthetic dentistry, teeth bleaching, low-level laser therapy (LLLT) and highintensity laser therapy (HILT), far more than the conventional lowpower diode laser.

- ✓ Low Speed Cutting
  - ✓ Coagulation
  - ✓ Hemostasis

#### 0.1-16W LASER

- ✓ Low Speed Cutting
  - ✓ Coagulation
  - ✓ Hemostasis
- ✓ High Speed Cutting
- ✓ Non-Charred Cutting
  - ✓ Teeth Bleaching
- ✓ High-Intensity Laser Therapy
- ✓ Low-Level Laser Therapy

# **SUPER-SHORT LASER PULSES REDUCE CHARRING** AND IMPROVE HEALING


reduce the charring of tissue and accelerate the postoperative healing. SOLASE Laser can generate a precise laser pulse - up to 10µs of pulse duration and 50kHz of



pulse duration=10µs



Continue wave mode

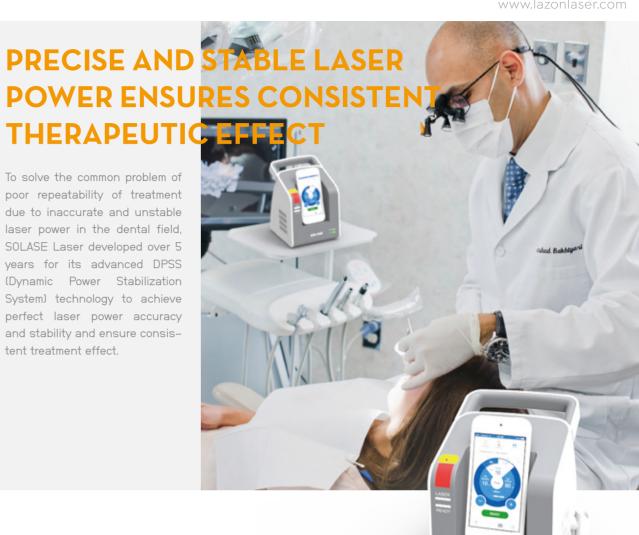

**EASY OPERATION:** AS SIMPLE AS A MASTER

#### **FOR BEGINNER**

- 30+ factory-loaded presets
- Step by step operating guide
- Easy installation and intuitive operation
- Plenty of reference material: video, literature, quick start, user manual...

#### **FOR PROFESSIONAL:**

- All parameters adjustable
- Patient records
- Pulse Interval or frequency display, Time or Energy display
- Time count-up or count-down

SOLASE Dental Diode Laser www.lazonlaser.com












## WHITENING HANDPIECE

WORLD'S FIRST FULL-MOUTH WHITENING HANDPIECE

Long and non-uniform laser irradiation will seriously increase pulp chamber temperature and may cause irreversible pulpal damage. LAZON provides the world's first full-mouth whitening handpiece to reduce the irradiation time to 1/4 of the conventional quarter mouth handpiece, with excellent uniform illumination to ensure same whitening effect on each tooth and preventing pulpal damage due to local intense illumination.

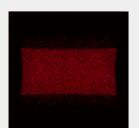

Laser spot of whitening handpiece from LAZON:
Large and Uniform




Laser spot of whitening handpiece from others:
Small and Non-Uniform

# THERAPY HANDPIECE

LARGER LIGHT SPOT COVERS LARGER AFFECTED AREA


The diameter of our new designed therapy handpiece can be adjusted from  $\Phi\,10\text{--}30\text{mm}$  to cover the affected area with different size.

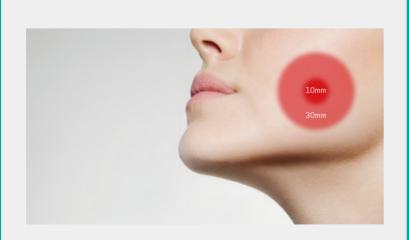

# **BIOSTIMULATION HANDPIECE**

DEEPER PENETRATION BY COLLIMATED LASER BEAM


Our specially designed biostimulation handpiece brings collimated laser beam, which can penetrate into the deeper tissue to get a better therapeutic effect.

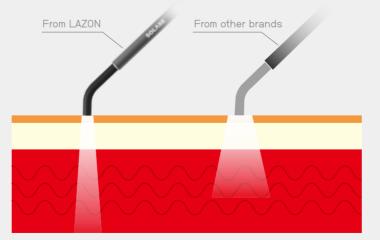

# **DISPOSABLE FIBER TIPS**

A NEW LENGTHENED DESIGN

We have creatively designed a lengthened disposable fiber tips that withstands autoclaving. It combines the tip and sleeve of traditional design to save your time and avoid cross-infection, without the needs of autoclaving the sleeve separately.

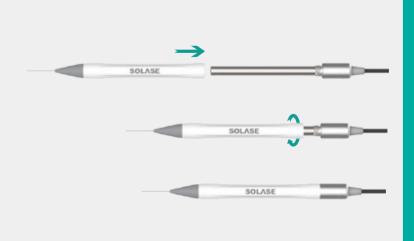

## **CUTTING FIBER HANDPIECE**

LOW COST AND DURABLE FROM BREAKING

Conventional reusable optical fiber has the inconvenience of fiber replacement, and the problem of being easily broken during use. We innovate our unique cutting fiber handpiece with low price but durable and easy to replace.

